

BMC Mini Auto CPAP System M1 Mini



en.bmc-medical.com



BMC Medical Co., Ltd.

3/F, Building No.4, No.1 Xinxing Road, Wuqing District, Tianjin 301700, China
Tel: +86-22-8293 9881 Fax: +86-22-8293 9881 etx. 810 en.bmc-medical.com

Specifications are subject to change without notice.

ST-SC-20-02-10-23-001-Ver1.0

tiny as me, free as you




Rich usage scenarios

BMC M1 Mini is an Auto-adjusting CPAP (Continuous Positive Airway Pressure) device with rich usage scenarios, designed for the treatment of Snoring and Obstructive Sleep Apnea Hypopnea Syndrome (OSAHS).





A variety of masks with
Moisture Exchanger available

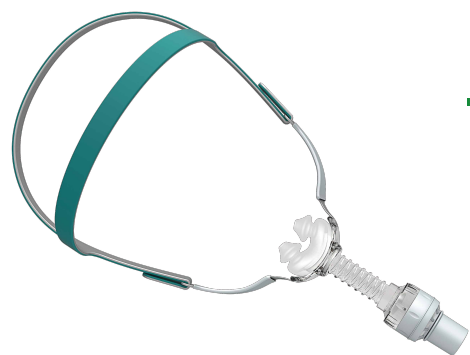

P2H Nasal Pillows Mask
Smaller contact area,
more comfortable experience.



N5H Nasal Mask
With light and smart groove
designed AutoFit Cushion,
it fits your face better.






N5AH Nasal Mask
Based on the forehead
frame-free design,
users can enjoy better vision.



BMC Mini App

M1 Mini adopts a minimalist style, with no screen design and touch buttons, showing high-end temperament. With M1 Mini App, you can easily—


-  Set up your device
-  Check the treatment effect
-  Upgrade the device firmware





Pressure adjusting sensitivity

Standard / Soft / Female,
three levels adjustable,
meet your personalized needs.

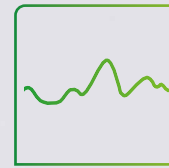

Smart A/C Feature

The intelligent pressure adjusting technology, the device can automatically adjust the treatment pressure scheme according to the changes of the patient's needs. It makes the treatment more suitable for you at any time.



Auto Ramp Feature

With the Auto ramp feature,
the device can sense
when the user falls asleep
and adjust the ramp time appropriately,
making the treatment
more comfortable and effective.



Accessories reminder

You can set the replacement reminders for the air filter,
moisture exchanger, mask and the air tubing,
so you can use M1 Mini without any worries.



Maintenance reminder

You can set a 180-day or 360-day maintenance reminder
to keep your M1 Mini in good condition.

Modern minimalist design, redefined “tiny”



Colorful ring light strip

Through different colors and flashes, it feeds back you a different experience.

Anti-slip design

There are four rubber anti-slip strips at the bottom, which effectively reduces the slip of the device during use.


Device size

Dimension: 159 mm x 66 mm x 72 mm

Weight: < 400 g



Product Accessories

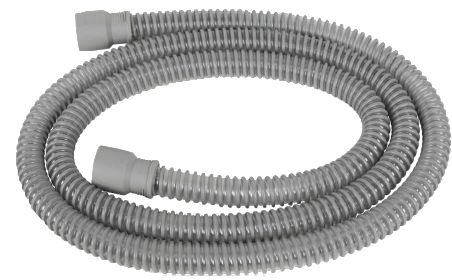

Air Tubing

- 15 / 22 mm optional



Carrying Case

- Reasonable space arrangement
- Convenient to carry
- Provide effective protection for M1 Mini



Storage Bag

- Soft material
- Compartment design
- Portable and convenient



Masks

- A variety of masks available
- Adapt to Moisture Exchanger



Air Filter

- Effectively filter dust and pollen in the air
- Please pay attention to regular cleaning and replacement



Moisture Exchanger

- Waterless humidification, keep your breath warm and moist
- 2 types optional, normal one and enhanced one
- Convenient and quick

Specifications

Scope of application	Suitable for the treatment of Snoring and Obstructive Sleep Apnea Hypopnea Syndrome (OSAHS).
Device Size	Dimensions: 159 mm × 66 mm × 72 mm Weight: < 400 g
Working Mode	CPAP, AutoCPAP
Pressure Range	4 - 20 hPa
Pressure Accuracy	± 0.5 hPa
Pressure Sensitivity	Standard / Soft / Female
Smart A/C	Yes
Ramp Time	0 - 60 min / Auto
Reslex (Expiration Pressure Release)	Off, Level 1 - 3
Auto On	Yes
Auto Off	Yes
Moisture Exchanger Type	Normal / Enhanced
Leak Alert	Yes
Accessories Reminder	Yes
Automatic Leakage Compensation	Yes
Firmware Upgrade through App	Yes
Wireless Data Transmission	Yes
Sound Pressure Level	< 30 dB, when the device is working at the pressure of 10 hPa.
AC Power Consumption	100 - 240 V AC, 50 / 60 Hz, 1.0 A max
Device Input	19 V, 1.26 A
Carrying Case	Optional

